

ACL7 SINGLE SHOT EPIDURAL - regional anesthesia kit
Clint Pharmaceuticals, Inc.

Drug Facts

APLICARE THREE POVIDONE-IODINE SWABSTICKS
[Aplicare, Inc.]

Three 4-inch saturated swabsticks

Povidone-iodine 10%

Antiseptic

Warnings

Do not use

- if allergic to iodine
- in the eyes

For external use only

Ask a doctor before use if injuries are

- deep or puncture wounds
- serious burns

Stop use and ask a doctor if

- redness, irritation, swelling or pain persists or increases
- infection occurs

Avoid pooling beneath patient

Keep out of reach of children. In case of accidental ingestion, seek professional assistance or consult a poison control center immediately.

NDC 52380-3101-5

NPN 02076101



← Tear Here Hold Upright Tear Here →

THREE

POVIDONE-IODINE SWABSTICKS

Three 4-inch saturated swabsticks

STERILE unless opened or damaged.

Drug Facts

Active ingredient	Purpose
Povidone-iodine USP 10%	Antiseptic

Use antiseptic skin preparation.

Warnings

Do not use if allergic to iodine or in the eyes

For external use only

Reorder No. S-3101

APLICARE EXP-01 2015

LOT

Drug Facts (continued)

Warnings

Ask a doctor before use if injuries are
■ deep or puncture wounds ■ serious burns

Stop use and ask a doctor if

■ redness, irritation, swelling or pain persists or increases
■ infection occurs

Avoid pooling beneath patient

Keep out of reach of children. In case of accidental ingestion, seek professional assistance or consult a poison control center immediately.

Directions apply locally as needed

Other information

■ 1% available iodine ■ latex free
■ for hospital or professional use only

Inactive ingredients citric acid, dicalcium phosphate, nonoxonyl-9, sodium hydroxide, water

For questions, comments, or to report serious side effects:

1-800-760-3236

Monday - Friday 8:30 a.m. - 5:00 p.m. EST



APLICARE, INC.
MERIDEN, CT 06450
U.S.A.
DOLLARD-des-ORMEAUX,
QUEBEC H9G 2S3
www.aplicare.com
1209

Package Label Display Panel



REF ACL7

To Reorder Call:
1-800-677-5022

PREP COMPONENTS:

- 1 Towel
- 1 Iodophor PVP Swabstick, 3's
- 3 Gauze Sponges

PROCEDURAL COMPONENTS:

- 1 Fenestrated Drape
- 1 Needle Stick Pad
- 1 18G X 1 1/2 in. Needle
- 1 25G X 1 1/2 in. Needle
- 1 Tuohy Epidural Needle Calibrated with Wings (20G X 3 1/2 in.)
- 1 Pulstar® L.O.R. Syringe (7ml, Luer Lock)
- 1 Plastic Syringe (10ml, Luer Lock)
- 1 Plastic Syringe (5ml, Luer Lock)
- 2 Plastic Syringes (3ml, Luer Lock)

SINGLE SHOT EPIDURAL

WARNINGS:

- A needle stick with a contaminated needle may cause infectious disease.
- The use of excessive force while placing needles into the stick pad may cause the needle to protrude through the bottom of the tray which may result in a contaminated needle stick.

PRECAUTIONS:

- **Use Aseptic technique.**
- To help prevent needle-stick injuries, needles should not be recapped or purposely bent. If excessive resistance is met during needle insertion, do not force the needle as damage may occur. To help avoid needle breakage, do not attempt to straighten a bent needle; discard it and complete the procedure with a replacement needle.
- After use, place sharps in a suitable sharps container. Dispose of contaminated product in a safe manner according to Centers for Disease Control and Prevention (USA) and Federal/State/Local regulations (EPA, OSHA) and health care facility guidelines or local equivalent.
- Do NOT Reuse/sterilize.
- To be used only by individuals familiar with single shot epidural procedures. For specific techniques and procedures, refer to standard textbooks.

Pulstar is a trademark of Smiths Medical. The symbol ® indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries.


Manufactured for:
Clint Pharmaceuticals
629 Shulte Lane
Old Hickory, TN 37138



Caution - Do Not Reuse - Does not contain Natural Rubber Latex - Do not use if package is damaged - Sterilized using ethylene oxide - Caution: Federal (U.S.A.) law restricts this device to sale by or on the order of a physician

Made in USA

PLTACL7 REV.003 11/11

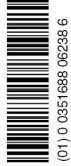

(01) 0 0351688 06238 6



REF ACL7

To Reorder Call:
1-800-677-5022

PREP COMPONENTS:

- 1 Towel
- 1 Iodophor PVP Swabstick, 3's
- 3 Gauze Sponges

PROCEDURAL COMPONENTS:

- 1 Fenestrated Drape
- 1 Needle Stick Pad
- 1 18G X 1 1/2 in. Needle
- 1 25G X 1 1/2 in. Needle
- 1 Tuohy Epidural Needle Calibrated with Wings (20G X 3 1/2 in.)
- 1 Pulstar® L.O.R. Syringe (7ml, Luer Lock)
- 1 Plastic Syringe (10ml, Luer Lock)
- 1 Plastic Syringe (5ml, Luer Lock)
- 2 Plastic Syringes (3ml, Luer Lock)

SINGLE SHOT EPIDURAL

WARNINGS:

- A needle stick with a contaminated needle may cause infectious disease.
- The use of excessive force while placing needles into the stick pad may cause the needle to protrude through the bottom of the tray which may result in a contaminated needle stick.

PRECAUTIONS:

- **Use Aseptic technique.**
- To help prevent needle-stick injuries, needles should not be recapped or purposely bent. If excessive resistance is met during needle insertion, do not force the needle as damage may occur. To help avoid needle breakage, do not attempt to straighten a bent needle; discard it and complete the procedure with a replacement needle.
- After use, place sharps in a suitable sharps container. Dispose of contaminated product in a safe manner according to Centers for Disease Control and Prevention (USA) and Federal/State/Local regulations (EPA, OSHA) and health care facility guidelines or local equivalent.
- Do NOT Reuse/sterilize.
- To be used only by individuals familiar with single shot epidural procedures. For specific techniques and procedures, refer to standard textbooks.

Pulstar is a trademark of Smiths Medical. The symbol ® indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries.

Manufactured for:
Clint Pharmaceuticals
629 Shulte Lane
Old Hickory, TN 37138



Caution - Do Not Reuse - Does not contain Natural Rubber Latex - Do not use if package is damaged - Sterilized using ethylene oxide - Caution: Federal (U.S.A.) law restricts this device to sale by or on the order of a physician

Made in USA

PLTACL7 REV.003 11/11



(01) 0 0351688 06238 6

ACL7 SINGLE SHOT EPIDURAL

regional anesthesia kit kit

Product Information

Product Type	MEDICAL DEVICE	Item Code (Source)	NHRIC:55553-454
---------------------	----------------	---------------------------	-----------------

Packaging

#	Item Code	Package Description	Marketing Start Date	Marketing End Date
1	NHRIC:55553-454-02	30 in 1 CASE		
1		1 in 1 PACKAGE, COMBINATION		

Quantity of Parts

Part #	Package Quantity	Total Product Quantity
Part 1	1 PACKET	6.5 mL

Part 1 of 1

APLICARE POVIDONE-IODINE TRIPLES

povidone-iodine solution

Product Information

Item Code (Source)	NDC:52380-3101
--------------------	----------------

Route of Administration	TOPICAL
-------------------------	---------

Active Ingredient/Active Moiety

Ingredient Name	Basis of Strength	Strength
POVIDONE-IODINE (UNII: 85H0HZU99M) (IODINE - UNII:9679TC07X4)	POVIDONE-IODINE	9.8 g in 96.4 mL

Inactive Ingredients

Ingredient Name	Strength
SODIUM PHOSPHATE, DIBASIC (UNII: GR686LBA74)	
CITRIC ACID MONOHYDRATE (UNII: 2968PHW8QP)	
SODIUM HYDROXIDE (UNII: 55X04QC32I)	
NONOXYNOL-9 (UNII: 48Q180SH9T)	
WATER (UNII: 059QF0K00R)	

Packaging

#	Item Code	Package Description	Marketing Start Date	Marketing End Date
1	NDC:52380-3101-5	6.5 mL in 1 PACKET		

Marketing Information

Marketing Category	Application Number or Monograph Citation	Marketing Start Date	Marketing End Date
OTC monograph not final	part333A	03/01/1998	

Marketing Information

Marketing Category	Application Number or Monograph Citation	Marketing Start Date	Marketing End Date
premarket notification	K965017	12/14/2005	

Labeler - Clint Pharmaceuticals, Inc. (609197785)**Registrant** - Smiths Medical ASD, Inc. (137835299)**Establishment**

Name	Address	ID/FEI	Business Operations
Smiths Medical ASD, Inc.		137835299	relabel, manufacture

Establishment

Name	Address	ID/FEI	Business Operations
------	---------	--------	---------------------

Aplicare, Inc.

107255002

manufacture

Revised: 6/2012

Clint Pharmaceuticals, Inc.